

## MODULAR STEEL LABORATORY CASEWORK

Modular steel casework is the best choice for any laboratory. It's not only long lasting durability, but also easy cleaning. Well maintained steel can serve decades without needing to be replaced. Our steel casework are with good design and perfect manufacturing; It is designed special for laboratory environments, chemical resistant, water resistant, scratch resistant.

Call now! +86-13302388256



## METAL LAB CASEWORK FEATURES

- Heavy duty, loading capacity 300kgs/m<sup>2</sup>;
- 100 lb. full extension drawer suspension;
- Famous brand hinges, soft and silence; 110 degree or 175 degree; more than 50,000 times, strong anti-rust and corrosion resistance, reaching the international hardware industry standard.


Phone: +86-13302388256

Email: [info@zhlabsfurniture.com](mailto:info@zhlabsfurniture.com)

Address: No.3068, Liangsha Road, Chendong Village, Zhongluotan Town, Baiyun District, Guangzhou City, Guangdong Province, China

- Modular design allows for simple field conversion from inset to overlay, door and drawer re-configurations and hardware component changes;
- Standard cabinet specification: W450mm-one door one drawer cabinet; W850mm-two doors two drawers cabinet;
- Each cabinet with drawers, doors and one adjustable shelf for easy storage;
- Doors and drawer panels are double wall structure, sound deadening materials are between panels;
- made entirely of 1.0mm thickness high quality cold rolled steel, welding structure, which provides strength and durability;
- Best paint treatments, to guarantee its durability.

**Standard Base Cabinet:** (Provided without locks)

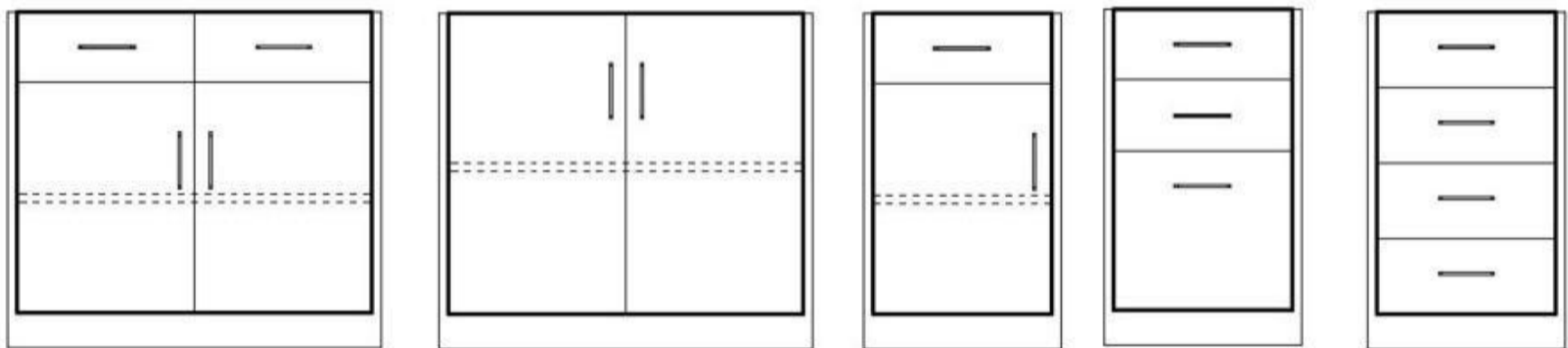

**Popular Color for Painted Steel Casework**



**PAINT TREATMENTS REQUIREMENTS**

- Pass through a phosphatizing process before painting.
- Coated with electrostatic application of epoxy resin powder and will be baked at high temperatures, to provide an uniform finish and resistance to chemicals as required in lab environments.

- The paint will be applied throughout the interior and exterior of the furniture. All components is individually painted, ensuring that there's no area vulnerable to corrosion due to lack of paint coverage.

At our local area, there are many raw material manufacturers and we keep good relationship with them. We have large and stable purchase quantity from them monthly. Cooperation with us, you can get high quality products with competitive price, and you will be always far away from clients complains. Contact our professional workers for assistant.

## STEEL LAB CASEWORK OPTIONAL

**Benchtop shelves:** The main frame is made of cold-rolled steel with 1.2mm thickness, coated with epoxy powder. The shelf is made of reinforcing glass with 12mm thickness and stainless steel tube to be guardrail. The shelves can be adjustable.

**Sink:** Stainless steel is popular in the international market; PP, epoxy resin, ceramic sink are available;

**Water faucet and gas tap;**

**Pegboard:** PP pegboard in black color, stainless steel pegboard in its color;

**Eyewash:** deck mounted eyewash beside water for first aid.

**Casters:** our lab workstations and tables can be moved and casters can be added.

Types	Accessories	Details	Loading Capacity	Material and Surface Finishing
Standard	Hinges	110 degree, soft closing (175 degree available)	100lbs	Stainless steel, stainless steel color
Standard	Slider	Not soft closing	100lbs	Zinc alloy, black or silver color
Special	Hinges	270 degree, not soft closing	100lbs	Zinc alloy, black color

Special	Slider	Heavy Bearing, Not soft closing	100lbs	Zinc alloy, black or silver color
---------	--------	------------------------------------	--------	--------------------------------------

## FAQ:

### What is the casework options?

Stainless steel, painted steel, wood, phenolic resin and specialty requests.

### Is the casework modular?

Yes, the modular design allows quickly shipment and easy maintenance in case there is any damage.

### What is the lab workstations and tables countertop made of?

Lab workstations and tables are available in different types of chemical-resistant tabletop materials such as epoxy resin, phenolic resin, stainless steel, and more.

### Are the lab workstations and tables movable?

Yes, our lab workstations and tables can be moved and casters can be added.

### Do I need to hire an architect and contractor?

No. We can provide the full design, layout, and build that an architect and contractor can offer. That saves you a lot of money!

### What's included in ZH LAB's design services?

Phone: +86-13302388256

Email: [info@zhlabsfurniture.com](mailto:info@zhlabsfurniture.com)

Address: No.3068, Liangsha Road, Chendong Village, Zhongluotan Town, Baiyun District, Guangzhou City, Guangdong Province, China



**ZH LABS FURNITURE**  
provide one-stop solution for your labs

<https://www.zhlabsfurniture.com/>

Our highly trained laboratory specialists can provide full design, CAD drawings, rivet files, architectural drawings, specifications, project estimates.

### **Do you offer OEM/ODM service?**

Yes, OEM/ODM service is available. Our experienced workers will fully understand the requirements and do the designs until you satisfied.

### **Which payment methods do you accept?**

T/T, Western Union. We will only ask for little part of deposit for the production arrangement, then we will offer pictures or video for quality checking. After the products are well packed and marked users name, the balance is needed. And you are welcome to visit us for quality checking.

**CONTACT US TODAY**

Phone: +86-13302388256

Email: [info@zhlabsfurniture.com](mailto:info@zhlabsfurniture.com)

Address: No.3068, Liangsha Road, Chendong Village, Zhongluotan Town, Baiyun District, Guangzhou City, Guangdong Province, China